

eMatrix™  
Complete  
Sublative  
Rejuvenation™



Science. Results. Trust.



Syneron®



"In my trials, I have found that the Sublative Rejuvenation™ treatment from Syneron offers a unique combination of ablation and skin resurfacing leading to wrinkle reduction and even lifting with less downtime than traditional skin rejuvenating treatments. My patients are pleased with the fast, visible results of smoother and brighter skin, without having to sacrifice time from their busy schedules."

Amy Taub, M.D.  
Dermatologist, Illinois

eMatrix™

Complete Sublative Rejuvenation™

The eMatrix system with Syneron's unique Sublative Rejuvenation treatment goes beyond fractional methods of skin resurfacing. The system utilizes Sublative RF™ – fractionated bi-polar radio frequency technology – to place the heat energy effectively into the dermis where it can produce significant dermal impact with minimal epidermal disruption.



# The Technology

*sub-la-tive re-ju-ve-na-tion*  
*[suh-bley-tiv ri-joo-vuh-ney-shuhn], noun.*

1. A new category of skin rejuvenation
2. High dermal impact, low epidermal disruption for all skin types
3. Minimal downtime, proven efficacy, satisfied patients
4. Practice-building procedure to generate life-long patients

# Deliver Effective Results With Minimal Downtime

The promise of noticeable results with minimal downtime and lower procedure cost is spurring interest in fractional methods of skin rejuvenation. Patients are increasingly requesting these treatments to improve their appearance.

The eMatrix system enables your practice to address this growing demand with beautiful results.

## Rejuvenate Your Practice and Meet the Growing Demand

The eMatrix system was designed with your practice in mind. It's an affordable system with a low total cost

of ownership. Adding Sublative Rejuvenation to your practice will generate greater patient traffic and revenue – leading to a strong ROI.

## Proven Clinical Efficacy



The results have been proven in clinical trials conducted by leading aesthetic physicians – and will assure that patients keep coming back and refer others to your practice.

## eMatrix System with the Sublative Rejuvenation Treatment – A Smart Choice

- Intelligent Feedback System™ provides real-time information on the energy delivery of each pulse
- Minimal downtime so patients can quickly resume their normal activities
- Deep dermal impact leading to the appearance of smoother skin – and resulting in satisfied patients
- Safe and effective for all skin types
- Quick treatments: typically 10 to 20 minutes depending on area treated; short treatment protocol: 3 treatments, 4-6 weeks apart
- Lightweight, compact and portable table-top design makes it easy to use anywhere in your practice – or transport between practice locations
- Excellent ROI potential – lower acquisition and disposable costs, and lowest overall cost of ownership compared to competitive products


## Complete Sublative Rejuvenation – High Density and Standard Sublative Rejuvenation Treatment

Dramatically improve texture, tone and wrinkles. Now offering two different low-cost, single-treatment disposable treatment tips – the original Sublative RF tip with 64 pins and the new Sublative RF High Density tip with 144 pins.

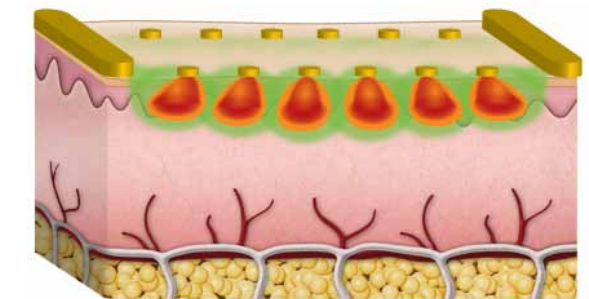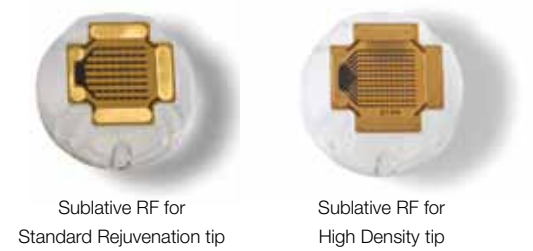

Sublative RF High Density tip with 144 pins

Photos: Tess Mauricio, M.D.

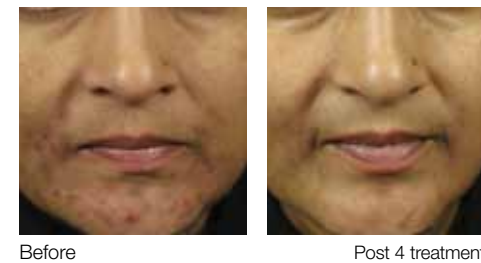

Photos: Sasha Jazayeri, M.D.




## Sublative Rejuvenation – Precise, Consistent Results

The eMatrix system's ergonomic applicator delivers bi-polar radio frequency energy via a matrix of electrodes of the disposable applicator tip. The RF energy generates fractional deep dermal heating in the region of the electrode matrix to induce skin injury, thus eliciting a wound healing response. Less affected tissue surrounding the matrix points may be stimulated by heat which helps accelerate the healing process, leading to a smoother, brighter appearance of the skin.


Only eMatrix features the Intelligent Feedback System to provide real-time feedback on the skin condition, impedance and exact energy delivered into the skin. This breakthrough technology helps assure consistency and reproducibility.

## SelectPulse™ – The Versatile Advantage

The SelectPulse feature of the eMatrix system gives clinicians the versatility to customize the Sublative Rejuvenation treatment for individual patients. You can select from three Sublative Rejuvenation programs to adjust the depth and density of dermal heating to treat a range of skin conditions.

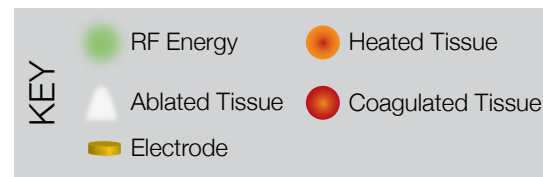

### Sublative Rejuvenation Program A

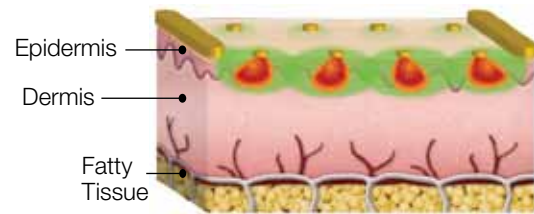

**Shallow Impact**  
Minor ablation, dermal heating and skin resurfacing – for improvements in overall skin tone and texture.

### Sublative Rejuvenation Program B




**Mid-Level Impact**  
Mild ablation, dermal heating and skin resurfacing – for color and textural irregularity correction.

### Sublative Rejuvenation Program C




**Deepest Impact**  
Moderate ablation and skin resurfacing.

## The Science of Sublative Rejuvenation

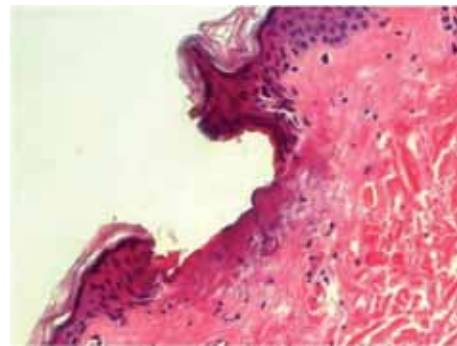

Histology showing immediate impact of Sublative RF Standard

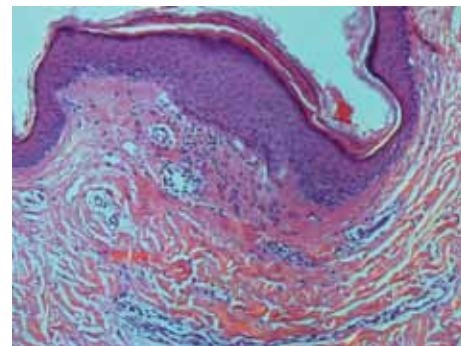

Histology showing healing process 3 days post treatment of Sublative RF Standard

## Sublative RF for the improvement of:

- Wrinkles
- Textural irregularities
- Superficial skin lesions





# Science. Results. Trust.



## System Specifications

Total RF energy	Up to 62 mJ/pin
Treatment programs	Program A, B or C
Applicator weight	0.5 kg / 1.1 lbs.
Applicator dimensions	15 x 16 cm / 5.9" x 6.3"
Disposable tip spot size	12 x 12 mm
Disposable Standard Sublative RF tip	64 pins
Disposable High Density Sublative RF tip	144 pins
Disposable treatment tip	Up to 200 pulses
User interface	15" touch screen
System dimensions	Table top: 53 x 30 x 40 cm / 21.8" x 11.8" x 15.7"
System weight	8.5 kg / 18.7 lbs.
Electrical requirements	100-230 VAC; 3A; 50-60 Hz



Portable table-top design.

[www.syneron.com](http://www.syneron.com) | [www.candelalaser.com](http://www.candelalaser.com)

Syneron Inc.  
3 Goodyear  
Suite A  
Irvine, California 92618  
USA  
Tel: 1 866 259 6661  
Email: [info@syneron.com](mailto:info@syneron.com)

Syneron Medical Ltd.  
Industrial Zone,  
Tavor Building  
P.O.B. 550 Yokneam Illit  
Israel 20692  
Tel: +972 73 244 2200  
Email: [info@syneron.com](mailto:info@syneron.com)

Syneron Medical (HK) Ltd.  
Room 3103, Hopewell Centre  
183 Queen's Road East  
Wanchai  
Hong Kong  
Tel: +852 25434326  
Email: [infoasia@syneron.com](mailto:infoasia@syneron.com)

Candela Corporation  
530 Boston Post Road  
Wayland, Massachusetts 01778  
USA  
Tel: 508 358 7400  
1 800 733 8550  
Email: [info@candelalaser.com](mailto:info@candelalaser.com)

Candela Iberica, SA  
Avda. de Castilla, 2  
Edificio Europa  
28830 Madrid  
Spain  
Tel: +34 91 656 85 63  
Email: [candela-iberica@candela-iberica.com](mailto:candela-iberica@candela-iberica.com)

